

MULTI-UNIT DASSYSTEM


Work different, work better



DYNAMIC ABUTMENT® SOLUTIONS


MULTI-UNIT DAS SYSTEM

Work different, work better

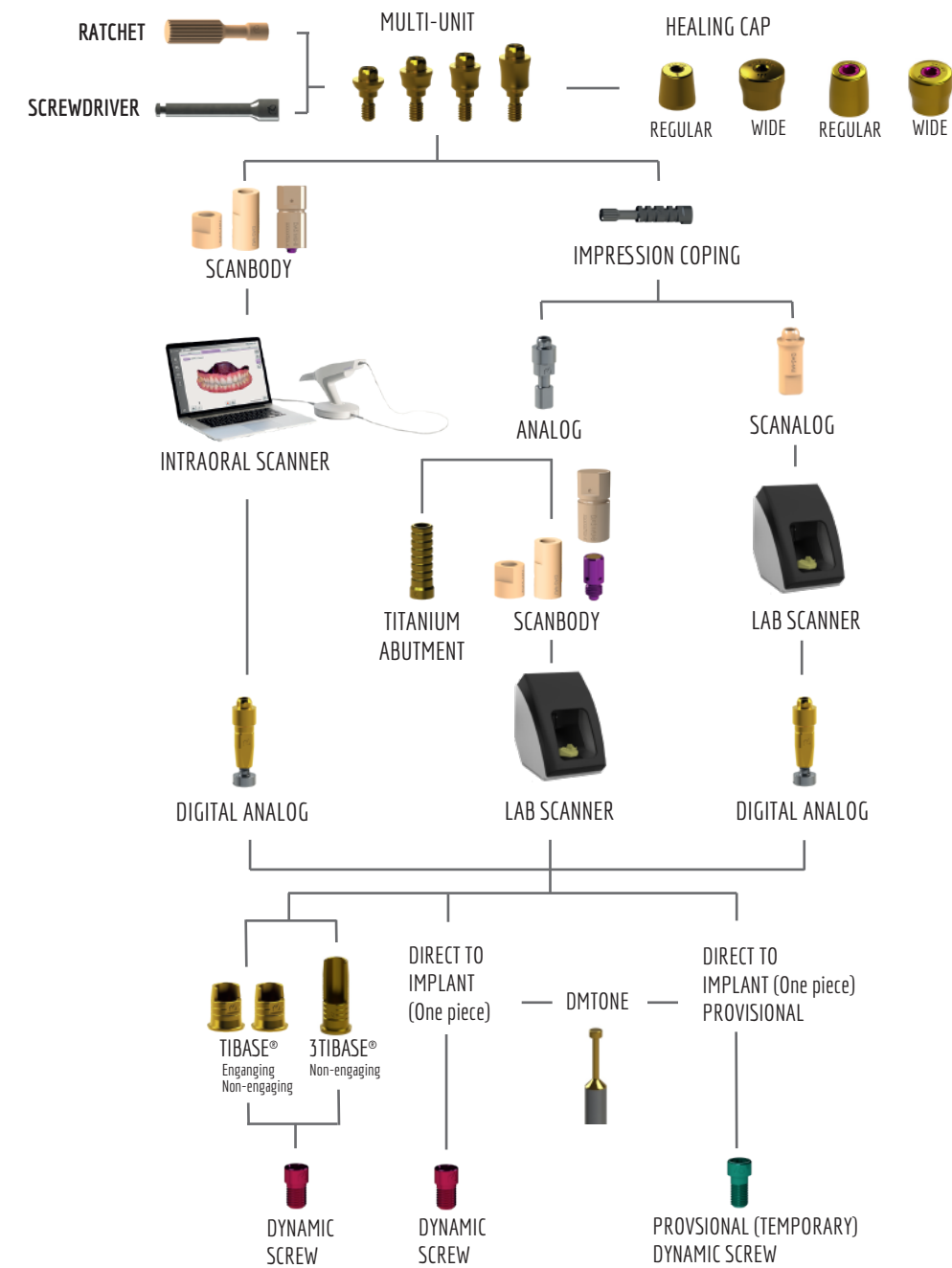

Digital workflow from the beginning to the end of the work. With both intraoral and lab scanners. Also available with full analogue workflow.


DIGITAL WORKFLOW - MULTI-UNIT DAS SYSTEM

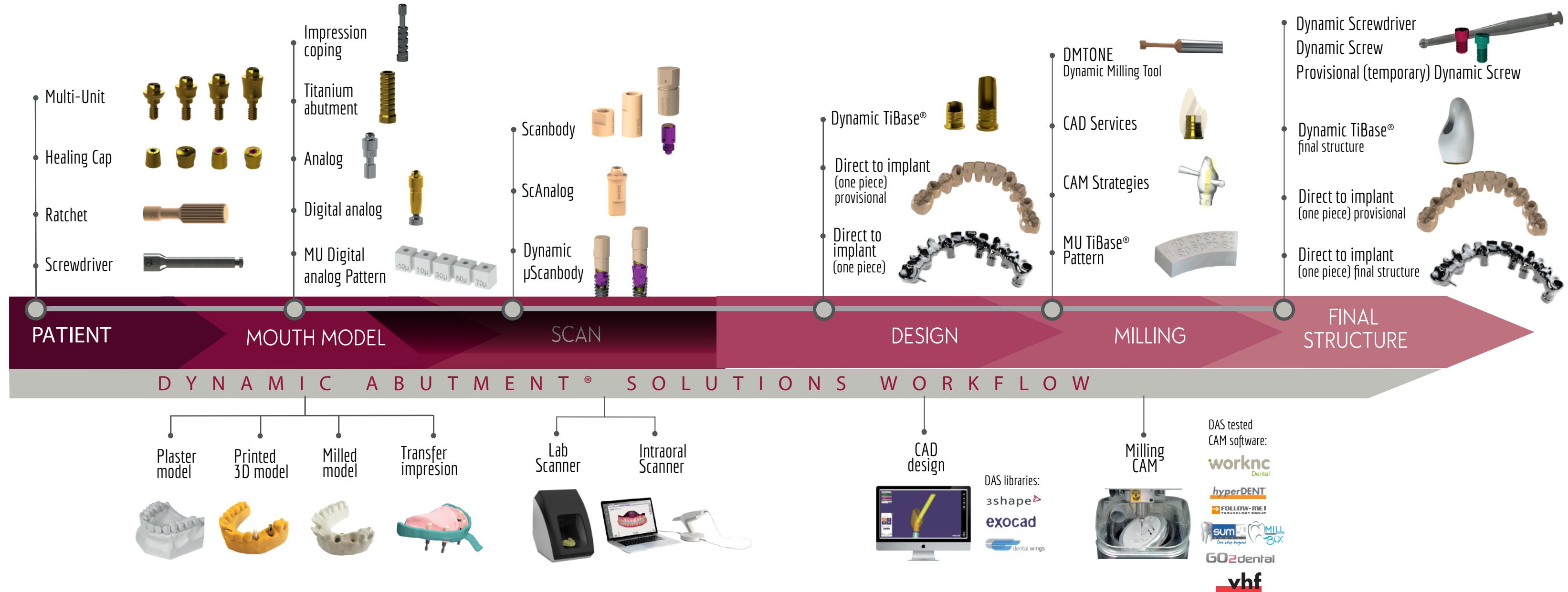

MULTI-UNIT

The Multi-unit abutment has been carefully designed to rehabilitate partially or fully edentulous arches, as well as individual.



Maximum diameter of the MU is 4mm, being narrower it avoids contact with the bone

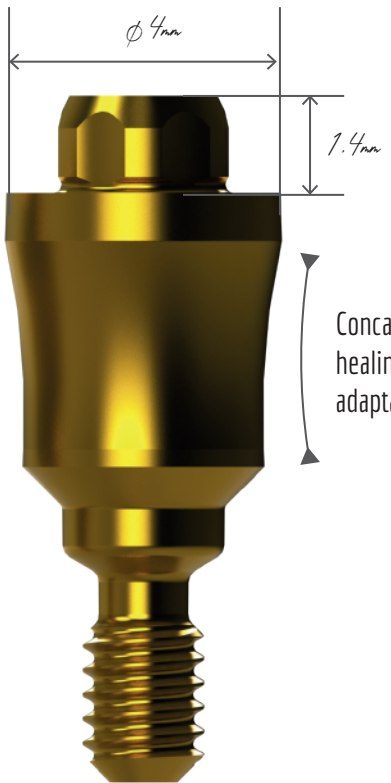


Lower cone height which allows greater disparelism than the classic multi-unit.
Divergence between implants: 50°

Concave design facilitates healing and soft tissue adaptation.


Available in 4 gingival heights.

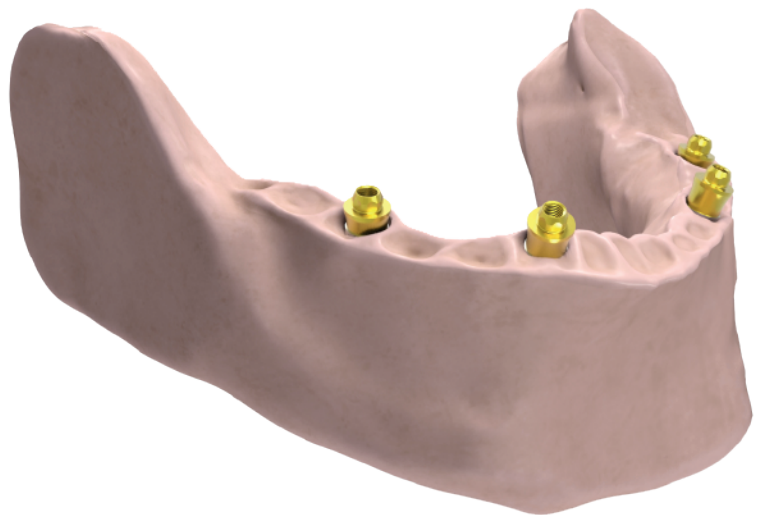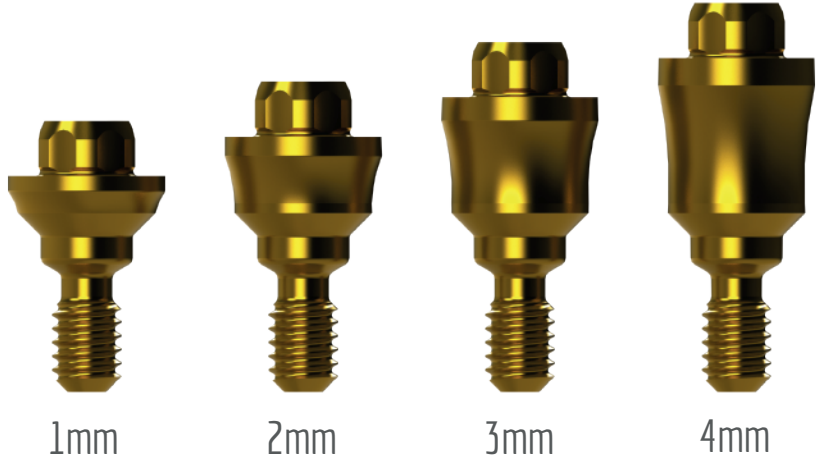


MU Insertion Procedure

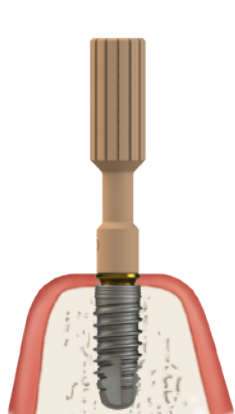

Step1
Use the plastic gripping tool to attach the Multi-Unit to the implant.


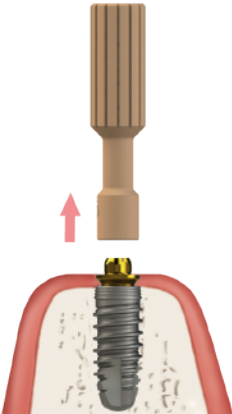


Step2
Remove the plastic tool.

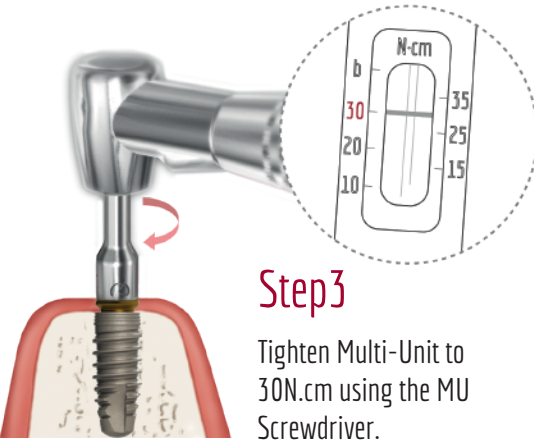


Step3
Tighten Multi-Unit to 30N.cm using the MU Screwdriver.



MU Healing Caps

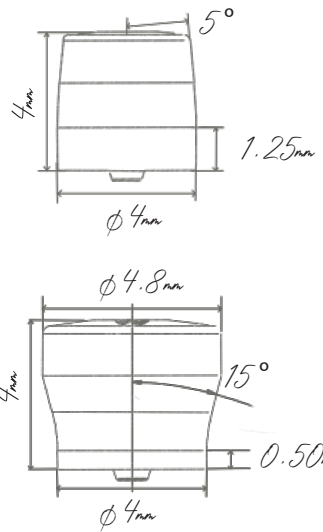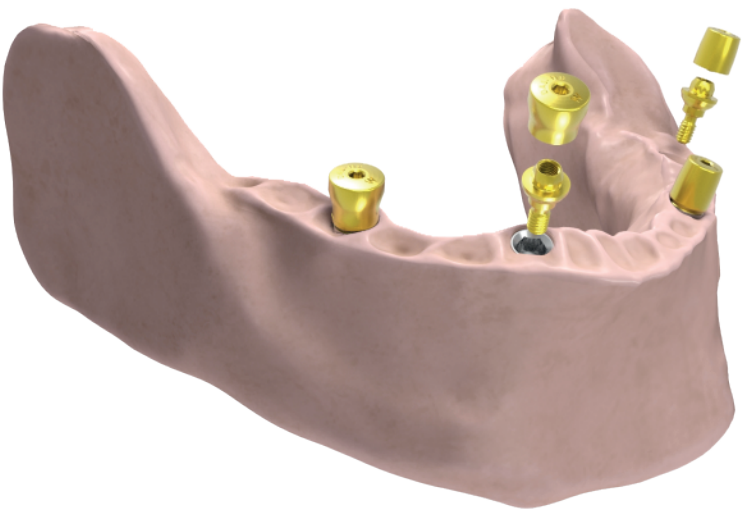

MULTI-UNIT SCANBODY

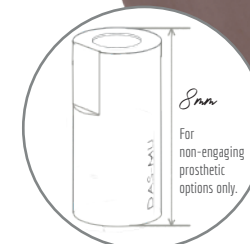


3 different scanbodies available in order to choose in which situation is better to use 4,5mm, 8mm or Dynamic μ Scanbody with magnet system.

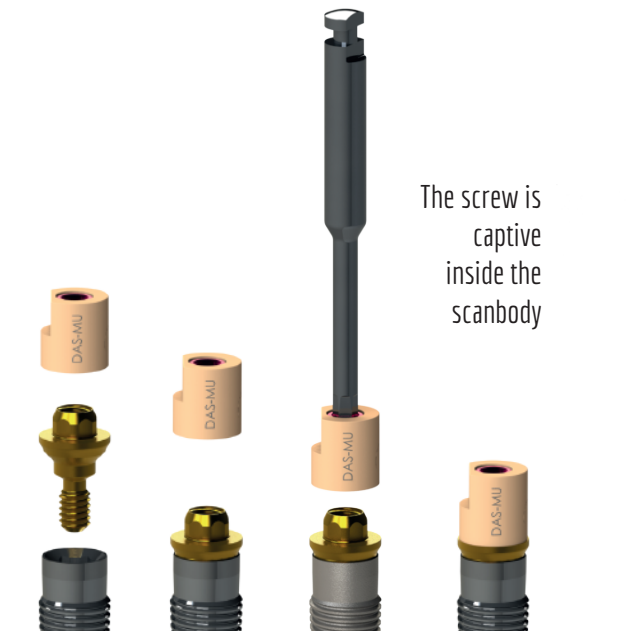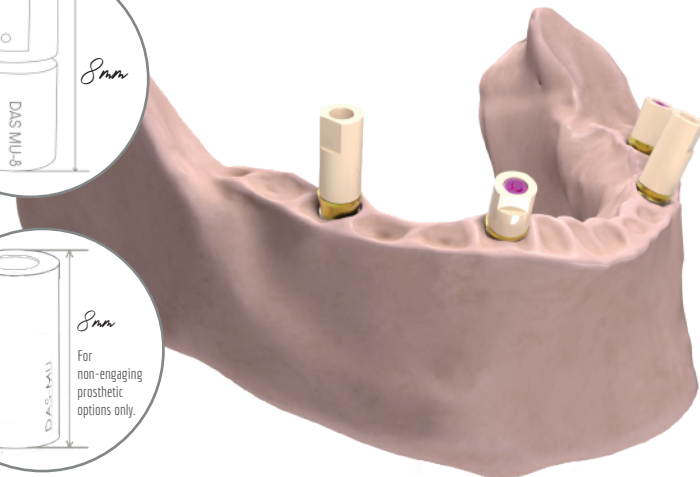

The screw is captive inside the scanbody

MULTI-UNIT SCANALOG

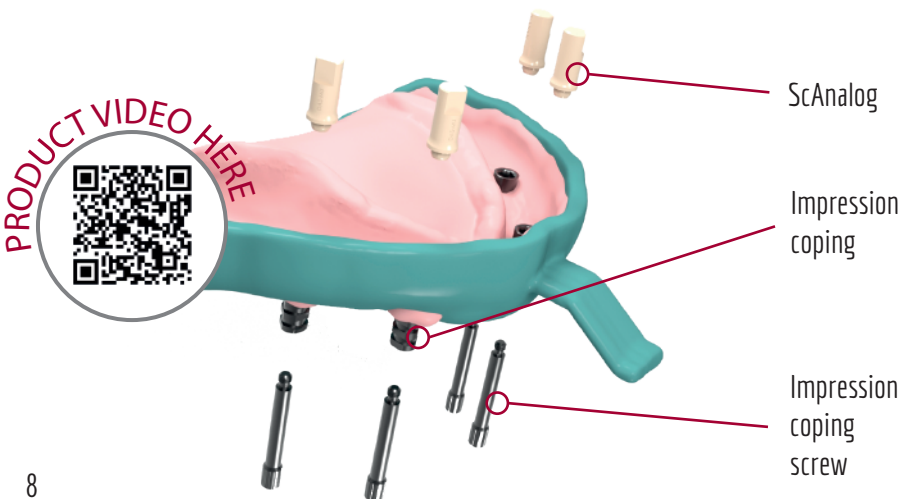


DYNAMIC μ SCANBODY

All components of the Multi-Unit DAS System can also be used with the Dynamic μ Scanbody to scan direct to implant and select virtually the ideal gingival height Multi-unit.


You can scan with the Dynamic μ Scanbody and then use the special Dynamic μ Scanbody library with the MU components.

For non-engaging prosthetic options only.

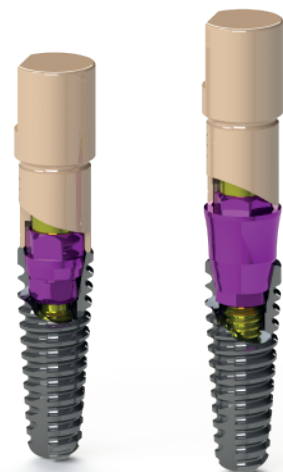

DIGITAL ANALOG


Digital analog of the dental implant to simulate implant position in a 3D printed dental model.

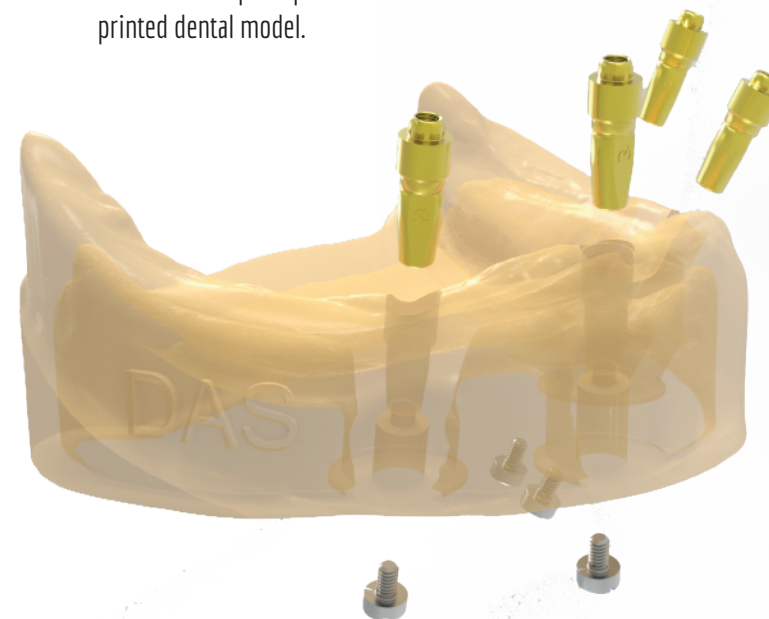

ANALOG



Also available traditional analog for plaster model.

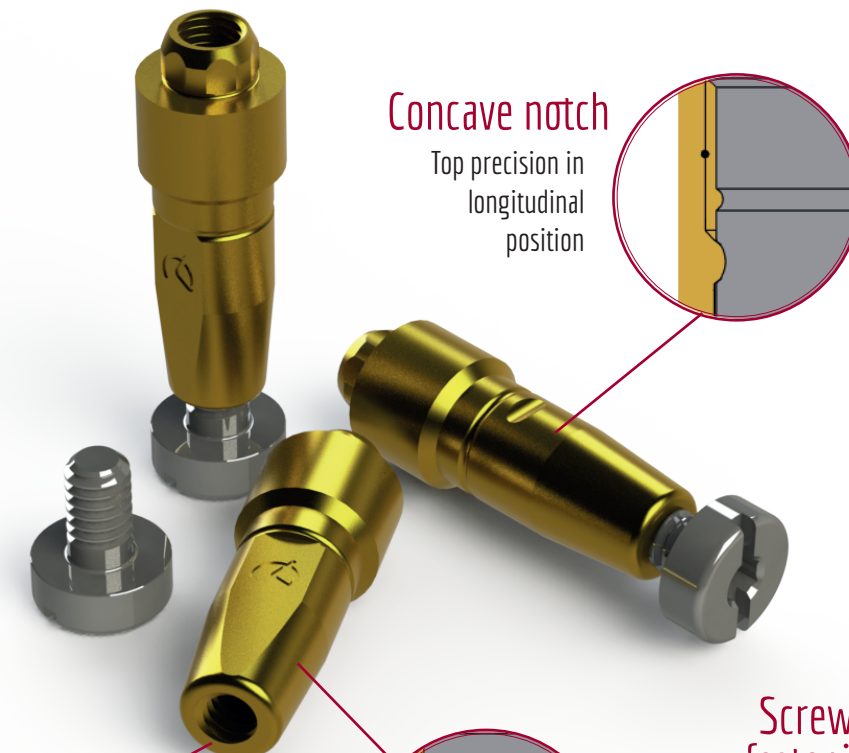

Concave notch
Top precision in longitudinal position

Screwed fastening
Prevents the analog from moving in Z

Curved Surface
Accuracy of orientation guaranteed


Longitudinal cut
Longitudinal cut to avoid rotation X-Y

DYNAMIC TIBASE®


To correct ANGULATION up to 45°




CAD DESIGN

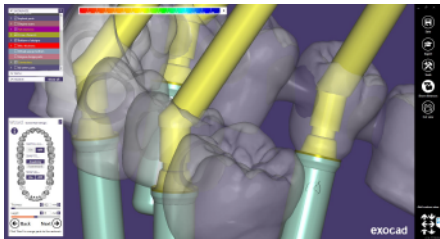

DAS Multi-U nit system **allows** to **change** the **MU** in the **library** without rescanning or redesigning the case. That facilitates the lab and clinic work, as technician can change the MU without the need to make a new appointment with the patient to re-scan.
If we make a NR case, the option of changing MU is not allow in the design.

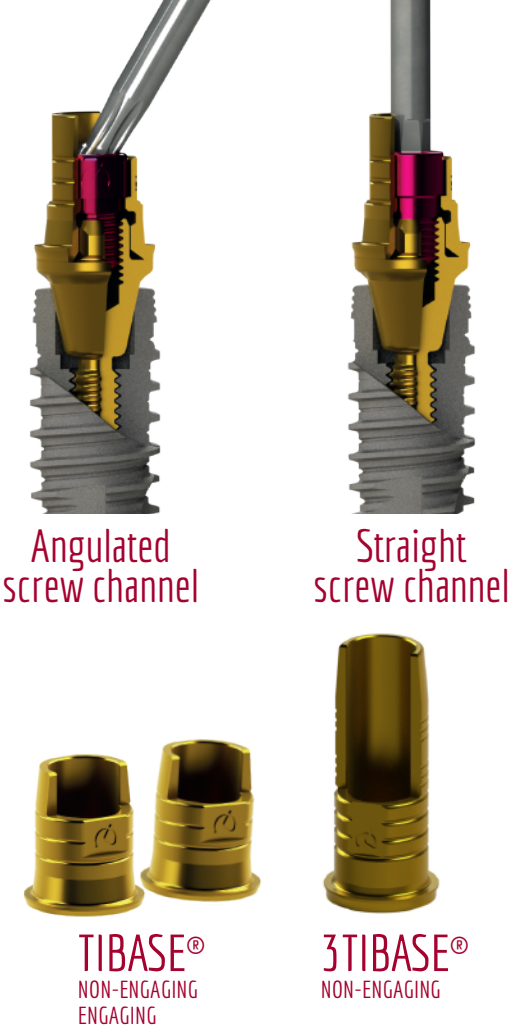


DIRECT TO IMPLANT (One piece)

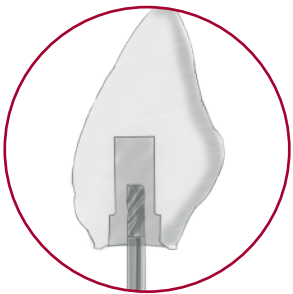

Step 1
Crown with pre drill



Step 2
Crown with angled chanel

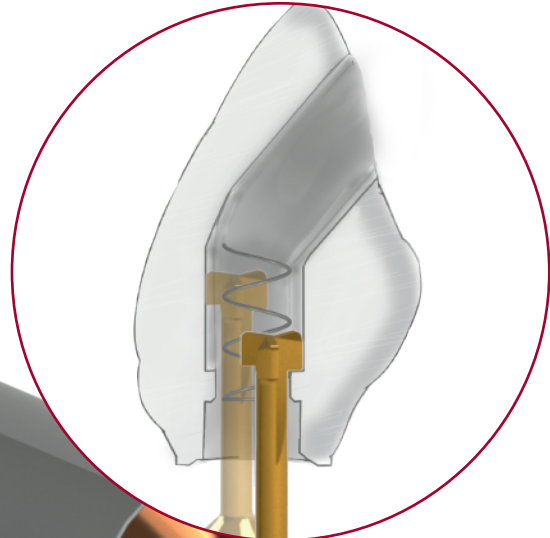

Step 3
DMTONE


Crown with Dynamic Milling Tool. Milling the screw seat and increasing the diameter of the straight channel.



Dynamic screw




Provisional (temporary) Dynamic screw

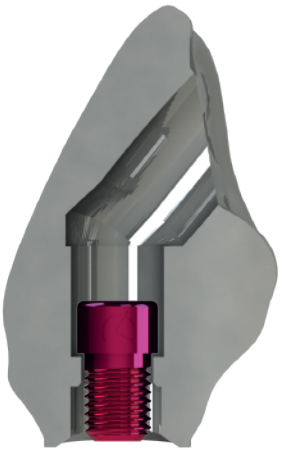

Final restorations 0-45°

Direct to implant libraries available R/NR

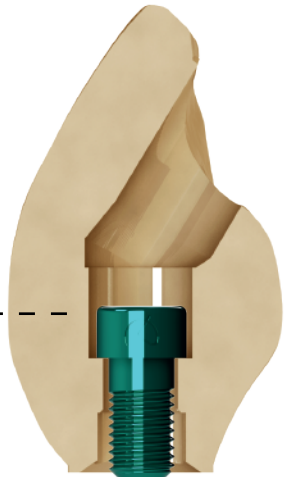


Provisional (temporary) restorations up to 35°

Provisional (temporary) Direct to implant libraries available R/NR

Greater seating thickness
Greater contact surface between the screw and the material for greater safety in soft materials.

DAS MU LIBRARIES

DAS MU PREMIUM library

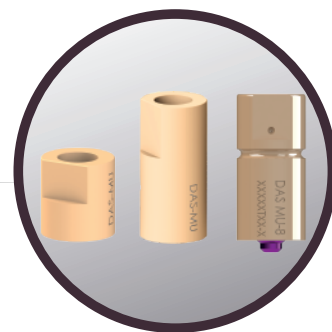


MU Scanbody

DAS_MU_4.5_XXXX*_PREMIUM
DAS_MU_8_XXXX*_PREMIUM
DAS_MU_DS8_XXXX*_PREMIUM



MU ScAnalog

DAS_MU_SA_XXXX*_PREMIUM

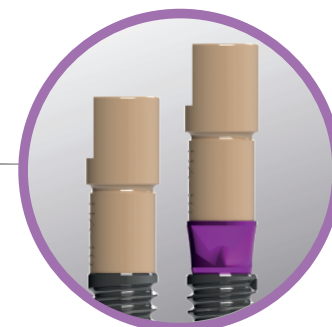

Dynamic μScanbody

DAS_MU_I/IG_XXXX*_PREMIUM

*XXXX= compatibility code

MU DAS SYSTEM COMPONENTS



Ratchet

49.409.000.01-2



Screwdriver

43.321.316.01-2



Healing Cap Regular

40.320.003.98-2



40.320.003.88-2 (Captive screw)



Healing Cap Wide

40.320.003.99-2



40.320.003.89-2 (Captive screw)



Impression coping

29.301.000.10-2 (Engaging)
29.301.000.11-2 (Non-engaging)



Analog

22.612.209.01-2



Titanium Abutment

35.312.209.21-2



Digital Analog

34.312.209.01-2



MU Scanbody 4,5mm

53.412.209.01-2



MU Scanbody 8mm

53.422.209.02-2 (R)



MU Dynamic μScanbody

52.418.137.01-2



Screwdriver adaptor

43.621.410.01-2
43.624.410.01-2
43.630.410.01-2



MU ScAnalog

23.412.209.01-2



MU Dynamic TiBase®

31.312.209.01-2 (Engaging)



31.322.209.01-2 (Non-engaging)



MU Dynamic 3TiBase®

31.322.209.21-2 (Non-engaging)



Dynamic Screw

41.320.040.01-2



Provisional (temporary) Dynamic Screw

41.320.050.02-2



Straight Screw

40.320.003.06-2



Dynamic Screwdriver

43.618.201.01-2 (18mm)
43.624.201.01-2 (24mm)
43.632.201.01-2 (32mm)




Screwdriver Hex.1,2

43.601.103.02-2



MU DMTone

33.390.716.01-2 Shank Ø3
33.490.716.01-2 Shank Ø4
33.690.716.01-2 Shank Ø6

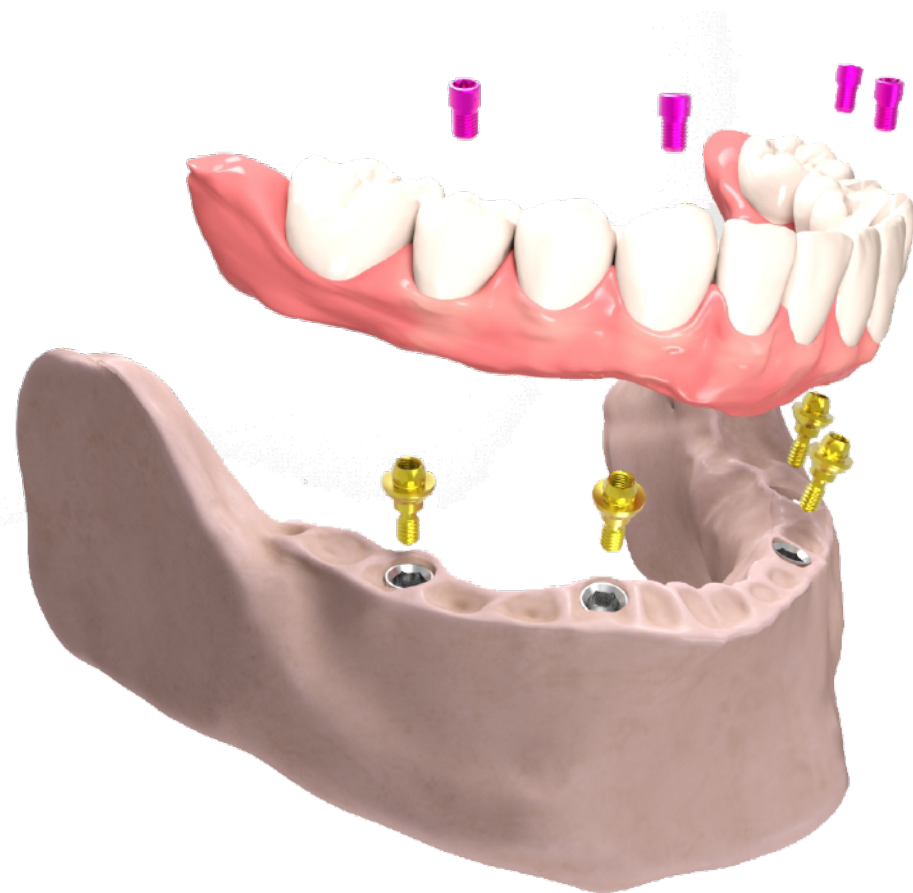

LIST OF COMPATIBILITIES AVAILABLE

AB	DENTIUM	NEOBIOTECH
ACE	DIO IMPLANTS	NOBEL BIO CARE
ADIN	EASY IMPLANT	NORIS MEDICAL
ALPHABIO	EUROTEKNIKA	OSTEOPLUS
ANCLADEN	HAHN IMPLANT (GLIDEWELL)	OSSTEM IMPLANT
ANKYLOS	HI-TEC	PALTOP
ARDS	HIOSEN	PROCLINIC
ASTRA	IBS	SEWON MEDIX
AVINENT	IDO IMPLANTS	SGS
BEGO	IHDE DENTAL (IMBIODENT)	SIGNO VINCES
BIOCONCEPT	IMPLANT DIRECT	SOUTHERN IMPLANTS
BIOGENESIS	IMPLANT GENESIS	STRAUMANN
BIOHORIZONS	INTRA-LOCK	TITANIUM - FIX
BIOMET 3i	JDENTALCARE	TRE-OSS
BIONER	LEADER	TRI DENTAL IMPLANTS
BIOTEC	MEDENTIS	UFIT
BTK	MEGAGEN	VULKAN IMPLANTS
B&W	MIS	YES IMPLANT
CAMLOG	MOZO-GRAU	ZIMMER
CONEXÃO SISTEMA DE PRÓTESE	MPI	
CORTEX	NEODENT	

AB

✱ I2	Implant: Ø 3,5/3,75/4,2/4,5/ 5/6 Platform: Standard (Code 0040)	p.26
✱ I22	Implant: Ø 3,75/4,22 Platform: Standard (Code 0040)	p. 26
✱ I5	Implant: Ø 3,5/3,75/4,2/4,5/5/6/7/8 Platform: Standard (Code 0040)	p. 26
✱ I55	Implant: Ø 3,75/4,2/4,5/5/6/7/8 Platform: Standard (Code 0040)	p. 26
✱ I10	Implant: Ø 4,2/5 Platform: Standard (Code 0040)	p. 26
✱ I15	Implant: Ø 6/7/8 Platform: Standard (Code 0040)	p. 26

ACE

✱ Infinity Internal Hex	Implant: Ø 3,7/4,1 Platform: 3,5 (Code 0040)	p. 26
-------------------------	---	-------

ADIN

✱ Swell	Implant: Ø 3,3 Platform: 3,45 (Code 0040)	p. 26
	Implant: Ø 3,75/4,2 Platform: 3,6 (Code 0040)	p. 26
	Implant: Ø 5 Platform: 4 (Code 0040)	p. 26
	Implant: Ø 6 Platform: 4,6 (Code 0040)	p. 26

✱ Touareg-S / Touareg-OS

	Implant: Ø 3,5 Platform: 3,45 (Code 0040)	p. 26
	Implant: Ø 3,75/4,2 Platform: 3,6 (Code 0040)	p. 26
	Implant: Ø 5 Platform: 4 (Code 0040)	p. 26
	Implant: Ø 6 Platform: 5 (Code 0040)	p. 26
✱ Touareg CloseFit	Implant: Ø 3,5 Platform: RP (Code 0021)	p. 24

ALPHABIO

✱ Internal Hex Connection (IH) SPI	Implant: Ø 3,3/3,75/4,2/5/6 Platform: Universal (Code 0040)	p. 26
✱ Internal Hex Connection (IH) ICE	Implant: Ø 3,7/3,75/4,2/4,65/5,3 Platform: Universal (Code 0040)	p. 26
✱ Internal Hex Connection (IH) DFI	Implant: Ø 3,3/3,75/4,2/4,5 Platform: Universal (Code 0040)	p. 26
✱ Internal Hex Connection (IH) ATID	Implant: Ø 3,3/3,75/4,2/5/6 Platform: Universal (Code 0040)	p. 26
✱ Internal Hex Connection (IH) NEO	Implant: Ø 3,75/4,2/5 Platform: 3,5 (Code 0040)	p. 26
✱ Conical Standard Connection (CS)	Implant: Ø 3,75/4,2/5 Platform: Standard (Code 0169)	p. 29

ANCLADEN

✱ Anclalock	Implant: Ø 3,75/4,25/5 Platform: 3,5 (Code 0040)	p. 26
-------------	---	-------

ARDS

✱ Smart	Implant: Ø 3,75/4,2/4,5 Platform: 3,75 (Code 0040)	p. 26
✱ Classic	Implant: Ø 3,75/4,2/4,5 Platform: 3,75 (Code 0040)	p. 26
	Implant: Ø 3,3/3,75/4,2/5/6 Platform: 3,75 (Code 0040)	p. 25
✱ Premium	Implant: Ø 3,3/3,75/4,2/5/6 Platform: 3,75 (Code 0040)	p. 26
✱ CIT	Implant: Ø 3,3/3,75/4,2/5/6 Platform: 3,75 (Code 0040)	p. 26

ANKYLOS

✱ Ankylos	Implant: Ø 3,5 Platform: 3,5 (Code 0075)	p. 26
	Implant: Ø 4,5 Platform: 4,5 (Code 0075)	p. 26
	Implant: Ø 5,5 Platform: 5,5 (Code 0075)	p. 26
	Implant: Ø 7 Platform: 7 (Code 0075)	p. 26

ASTRA

✱ Aqua	Implant: Ø 3,5/4 Platform: Aqua (Code 0004)	p. 22
✱ Evolution (Internal)	Implant: Ø 3,6 Platform: 3,6 (Code 0006)	p. 23
	Implant: Ø 4,2 Platform: 4,2 (Code 00007)	p. 23
	Implant: Ø 4,8 Platform: 4,8 (Code 00091)	p. 27

AVINENT

✱ HI/IC	Implant: Ø 3,1//3,5/4 Platform: 3,5 (Code 0040)	p. 26
---------	--	-------

BIOCONCEPT

✱ BC Bone Level	Implant: Ø 3,3 Platform: Narrow (Code 0033)	p. 25
	Implant: Ø 4,1/4,8 Platform: Regular (Code 0035)	p. 26
✱ BV Tapered Bone Level	Implant: Ø 3,5 Platform: Narrow (Code 0029)	p. 25
	Implant: Ø 4/4,5/5 Platform: Regular (Code 0030)	p. 25

BIOHORIZONS

✱ Tapered Internal	Implant: Ø 3,8 Platform: 3,5 (Yellow) (Code 0040)	p. 26
	Implant: Ø 3/3,4 Platform: 3,5 (Grey) (Code 0102)	p. 27
✱ Internal	Implant: Ø 3,5/4 Platform: 3,5 (Yellow) (Code 0040)	p. 26

BIOMET 3i

✱ Certain	Implant: Ø 3,25/4 Platform: 3,4 (Code 0001)	p. 22
	Implant: Ø 4/5 Platform: 4,1 (Code 0001)	p. 22

BIONER

✱ Hikelt	Implant: Ø 3,8 Platform: 3,95 (Code 0040)	p. 26
✱ TopDM	Implant: Ø 3,5 Platform: 3,5 (Code 0021)	p. 24
	Implant: Ø 4 Platform: 4 (Code 0021)	p. 24
	Implant: Ø 5 Platform: 5 (Code 0021)	p. 24

BIOTEC

✱ SPR/CIM	Implant: Ø 3,3 Platform: 3,3 (Code 0040)	p. 26
	Implant: Ø 3,75 Platform: 3,75 (Code 0040)	p. 26
✱ SPR/SPTT/CIM	Implant: Ø 4,2 Platform: 4,2 (Code 0040)	p. 26
	Implant: Ø 5 Platform: 5 (Code 0040)	p. 26

BTK

✱ Klassic / Konic	Implant: Ø 3,25/4 Platform: 3,5 IR (Code 0040)	p. 26
-------------------	---	-------

B&W

✱ Internal Hexagon CIH	Implant: Ø 3,3/4 Platform: 4 (Code 0040)	p. 26
------------------------	---	-------

CAMLOG

✱ Conelog Screw-Line	Implant: Ø 3,8 Platform: 3,8 (Code 0120)	p. 27
	Implant: Ø 4,3 Platform: 4,3 (Code 0121)	p. 28

CONEXÃO SISTEMA DE PRÓTESE

✱ Flash	Implant: Ø 3,5/4,3/5 Platform: Universal (Code 0021)	p. 24
✱ Torq	Implant: Ø 3,5/3,75/4 Platform: Universal (Code 0021)	p. 24
✱ Expand	Implant: Ø 3,75/4/5 Platform: Universal (Code 0021)	p. 24

CORTEX

✱ Dynamix	Implant: Ø 3,3/3,8/4,2/5/6 Platform: 3,75 (Code 0040)	p. 26
✱ Classix	Implant: Ø 3,3/3,8/4,2/5/6 Platform: 3,75 (Code 0040)	p. 26
✱ Saturn	Implant: Ø 3,8/4,2 Platform: 3,5 (Code 0040)	p. 26
✱ Conical Platform:	Implant: Ø 3,3/3,8/4,2 Platform: RP (Code 0004)	p. 22
✱ Magix	Implant: Ø 3,3/3,8/4,2 Platform: RP (Code 0004)	p. 22

DENTIS

✿ s-Clean	
Implant: Ø 3,7	
Platform: Mini (Code 0030)	p. 25
Implant: Ø 4,1/4,3	
Platform: Regular (Code 0030)	p. 25
Implant: Ø 4,8	
Platform: Wide (Code 0030)	p. 25

DENTIUM

✿ SuperLine and Implantium	
Implant: Ø 3,4	
Platform: 3,6 (Code 0030)	p. 25
Implant: Ø 3,8	
Platform: 4 (Code 0030)	p. 25
Implant: Ø 4,3	
Platform: 4,5 (Code 0030)	p. 25
Implant: Ø 4,8	
Platform: 5 (Code 0030)	p. 25
Implant: Ø 4,8	
Platform: 6 (Code 0030)	p. 25

DIO IMPLANTS

✿ UF II	
Implant: Ø 3,8/4/4,5/5/5,5	
Platform: Regular (Code 0030)	p. 25
✿ UF II Narrow	
Implant: Ø 3 / 3,3	
Platform: Narrow (Code 0014)	p. 23

EASY IMPLANT

✿ Master C	
Implant: Ø 3,5	
Platform: 3,5 (Ocean) (Code 0004)	p. 22
Implant: Ø 4	
Platform: 4 (Ocean) (Code 0004)	p. 22

Implant: Ø 4,5	
Platform: 4,5 (Lilas) (Code 0030)	p. 25
Implant: Ø 5	
Platform: 5 (Lilas) (Code 0030)	p. 25

✿ Master S

Implant: Ø 3,3	
Platform: 3,3 (Ocean) (Code 0004)	p. 22
Implant: Ø 3,75	
Platform: 3,75 (Lilas) (Code 0030)	p. 25
Implant: Ø 4,25	
Platform: 4,25 (Lilas) (Code 0030)	p. 25
Implant: Ø 4,75	
Platform: 4,75 (Lilas) (Code 0030)	p. 25

✿ Master L

Implant: Ø 3,3	
Platform: 3,3 (Lilas) (Code 0030)	p. 25
Implant: Ø 3,75	
Platform: 3,75 (Lilas) (Code 0030)	p. 25
Implant: Ø 4,25	
Platform: 4,25 (Lilas) (Code 0030)	p. 25
Implant: Ø 4,75	
Platform: 4,75 (Lilas) (Code 0030)	p. 25

EUROTEKNIKA

✿ Naturactis	
Implant: Ø 3,5	
Platform: 3,4 (Code 0004)	p. 22
Implant: Ø 4	
Platform: 3,8 (Code 0004)	p. 22
Implant: Ø 4,5	
Platform: 4,3 (Code 0004)	p. 22
Implant: Ø 5	
Platform: 4,8 (Code 0004)	p. 22

✿ Natea

Implant: Ø 3,6/4,1/4,8	
Platform: Narrow (Code 0004)	p. 22
Implant: Ø 3,6/4,1/4,8	
Platform: Regular (Code 0004)	p. 22
Implant: Ø 6	
Platform: Wide (Code 0004)	p. 22

✿ Naturall

Implant: Ø 3,5	
Platform: Narrow (Code 0004)	p. 22
Implant: Ø 4/4,5	
Platform: Regular (Code 0004)	p. 22
Implant: Ø 5	
Platform: Wide (Code 0004)	p. 22

GALIMPLANT

✿ Internal Connection

Implant: Ø 3,5	
Platform: 3,5 (Code 0004)	p. 22
Implant: Ø 4	
Platform: 4 (Code 0004)	p. 22
Implant: Ø 5	
Platform: 5 (Code 0004)	p. 22

HAHN IMPLANT (GLIDEWELL)

✿ Hahn Tapered Implant	
Implant: Ø 3,5/4,3	
Platform: 3,5/4,3 (Code 0021)	p. 24
Implant: Ø 5	
Platform: 5 (Code 0022)	p. 24

HI-TEC

✿ Tapered Self Thread	
Implant: Ø 3,3/3,75	
Platform: 3,5 (Code 0040)	p. 26

✿ Logic Plus

Implant: Ø 3,5	
Platform: 3,7 (Code 0040)	p. 26
Implant: Ø 4,3	
Platform: 3,9 (Code 0040)	p. 26

HIOSSEN

✿ ETTI SA/ETIII SA

Implant: Ø 3,5	
Platform: Mini (Code 0029)	p. 25
Implant: Ø 4/4,5/5	
Platform: Regular (Code 0030)	p. 25

✿ ETTI BA

Implant: Ø 3,5	
Platform: Mini (Code 0029)	p. 25
Implant: Ø 4/4,5/5	
Platform: Regular (Code 0030)	p. 25

IBS

✿ Magic FC

Implant: Ø 4/4,5/5/5,5/6/6,5	
Platform: 3,8 (Code 0030)	p. 25

✿ N.R Fix

Implant: Ø 3/3,5	
Platform: 3,8 (Code 0030)	p. 25

IDO IMPLANTS

✿ IDo Implant

Implant: Ø 3,8/4/4,5/5/5,5/6/7	
Platform: Universal (Code 0030)	p. 25

IHDE DENTAL (IMBIODENT)

✿ Bone Level Plus

Implant: Ø 3,3	
Platform: 3,3 (Code 0033)	p. 25
Implant: Ø 4,1	
Platform: 4,1 (Code 0035)	p. 26
Implant: Ø 4,8	
Platform: 4,8 (Code 0035)	p. 26

IMPLANT DIRECT

✿ Legacy

Implant: Ø 3,7/4,2	
Platform: 3,5 (Code 0040)	p.26

✿ SwishActive

Implant: Ø 3,3	
Platform: 3 (Code 0021)	p. 24
Implant: Ø 4,1/4,8	
Platform: 3,4 (Code 0022)	p. 24

✿ Interactive

Implant: Ø 3,2/3,7	
Platform: 3 (Code 0021)	p. 24
Implant: Ø 4,3/5	
Platform: 3,4 (Code 0022)	p. 24

IMPLANT GENESIS

✿ Aktiv System

Implant: Ø 3,5/3,75/4,2/5	
Platform: Standard (Code 0040)	p. 26

INTRA-LOCK

✿ IntraHex

Implant: Ø 3,75/4	
Platform: 3,5 (Code 0040)	p. 26

JDENTALCARE

✿ JDEvolution/JDEvolution Plus

Implant: Ø 3,7	
Platform: 3,7 (Code 0040)	p. 26
Implant: Ø 4,3/5	
Platform: 4 (Code 0040)	p. 26
Implant: Ø 6	
Platform: 5 (Code 0040)	p. 26

✿ JD ICON

Implant: Ø 3,9	
Platform: 3,9 (Code 0022)	p. 24
Implant: Ø 4,3	
Platform: 4 (Code 0022)	p. 24
Implant: Ø 5	
Platform: 4,7 (Code 0022)	p. 24

LEADER

✿ Tixos Internal Hex

Implant: Ø 3,3	
Platform: 3,5 (Code 0040)	p. 26
Implant: Ø 3,75	
Platform: 4 (Code 0040)	p. 26

MEDENTIS

✿ ICX - Templant

Implant: Ø 3,75	
Platform: 3,75 (Code 0125)	p. 28
Implant: Ø4,1	
Platform: 4,1 (Code 0125)	p. 28
Implant: Ø 4,8	
Platform: 4,8 (Code 0125)	p. 28

MEGAGEN

- AnyRidge
 - Implant: Ø 3,5
 - Platform: Small (Code 0015) p. 24
- Implant: Ø 4/4,5
 - Platform: Regular (Code 0015) p. 24
- Implant: Ø 5/5,5
 - Platform: Wide (Code 0015) p. 24
- AnyOne Internal
 - Implant: Ø 3,5/4/4,5/5/6/7
 - Platform: General (Code 0030) p. 25
- Mini Narrow Ridge
 - Implant: Ø 3 / 3,4
 - Platform: Mini (Code 0014) p. 23

MIS

- Seven
 - Implant: Ø 3,75/4,2
 - Platform: Standard (Code 0040) p. 26
- M4
 - Implant: Ø 3,75/4,2
 - Platform: Standard (Code 0040) p. 26

MOZO-GRAU

- MG Inhex
 - Implant: Ø 3,75/4,25
 - Platform: 2,8 Standard (Code 0004) p. 22

MPI

- Privilege CM
 - Implant: Ø 3,5/4
 - Platform: Regular (Code 0004) p. 22
- Excellence CM
 - Implant: Ø 3,5/4
 - Platform: Regular (Code 0004) p. 22

NEODENT

- Helix GM/Drive GM/ Titamax GM
 - Implant: Ø 3,5/3,75/4/4,3/5/6
 - Platform: Regular (Code 0186) p. 28

NEOBIOTECH

- IS Implant: System
 - Implant: Ø 4
 - Platform: Regular 4 (Code 0030) p. 25
- Implant: Ø 4,5
 - Platform: Regular 4,5 (Code 0030) p. 25
- Implant: Ø 5
 - Platform: Wide 5 (Code 0030) p. 25

NOBEL BIOCARE

- Active
 - Implant: Ø 3,5
 - Platform: Narrow (Code 0021) p. 24
- Implant: Ø 4,3/5
 - Platform: Regular (Code 0022) p. 24
- NobelParallel
 - Implant: Ø 3,75
 - Platform: Narrow (Code 0021) p. 24
- Implant: Ø 4,3/5
 - Platform: Regular (Code 0022) p. 24

NORIS MEDICAL

- Tuff
 - Implant: Ø 3,3/3,75/4,2/5/6
 - Platform: 3,75 (Code 0040) p. 25
- Tuff TT
 - Implant: Ø 3,3/3,75/4,2/5/6
 - Platform: 3,75 (Code 0040) p. 26
- Onix
 - Implant: Ø 3,3/3,75/4,2/5/6
 - Platform: 3,75 (Code 0040) p. 26

- Cortical
 - Implant: Ø 4,0/5/6
 - Platform: 3,75 (Code 0040) p. 26

- PteryCore
 - Implant: Ø 4,2
 - Platform: 3,75 (Code 0040) p. 26

- PteryFit
 - Implant: Ø 4,2
 - Platform: 3,75 (Code 0040) p. 26

OSTEOPLUS

- Shi
 - Implant: Ø 3,3 / 3,75 / 4,2
 - Platform: 3,5 (Code 0040) p. 26

OSSTEM IMPLANT

- TS
 - Implant: Ø 3,5
 - Platform: Mini (Code 0029) p. 25
- Implant: Ø 4/4,5/5/6/7
 - Platform: Regular (Code 0030) p. 25

PALTOP

- Conical Active
 - Implant: Ø 3,75/4,2/5
 - Platform: Standard (Code 0029) p. 25

PROCLINIC

- Aqua CM
 - Implant: Ø 3,5/4/5
 - Platform: 2,82 (Code 0004) p. 22
- Cylindrical Internal/Conical Internal
 - Implant: Ø 3,3/3,75/4,25/5//3,5/4/5
 - Platform: 3,5 (Code 0040) p. 26

SEWON MEDIX

- IH2 SLA SYSTEM
 - Implant: Ø 3,5
 - Platform: Mini (Code 0029) p. 25
- Implant: Ø 3,5/4/4,5/5
 - Platform: Regular (Code 0030) p. 25

- IH2 RBM SYSTEM
 - Implant: Ø 3,5
 - Platform: Mini (Code 0029) p. 25
- Implant: Ø 3,5/4/4,5/5
 - Platform: Regular (Code 0030) p. 25

SGS

- P1
 - Implant: Ø 3,2/3,75/4,2/5/6
 - Platform: 3,75 (Code 0040) p. 26
- P7
 - Implant: Ø 3,2/3,75/4,2/4,5/5/6
 - Platform: 3,75 (Code 0040) p. 26

SIGNO VINCES

- Compact
 - Implant: Ø 4,5
 - Platform: CM3,8 (Code 0004) p. 22
- Duocon
 - Implant: Ø 3,8
 - Platform: CM3,8 (Code 0004) p. 22
- Infra
 - Implant: Ø 3,3/3,8/4,6
 - Platform: CM (Code 0004) p. 22

SOUTHERN IMPLANTS

- Deep Conical
 - Implant: Ø 3,5/4
 - Platform: 2,95/3,1 (Code 0004) p. 22
- Internal Hex
 - Implant: Ø 3,75/4,2/5
 - Platform: Universal (Code 0040) p. 26

STRAUMANN

- Bone Level
 - Implant: Ø 3,3
 - Platform: NC- 3,3 (Code 0033) p. 25
- Implant: Ø 4,1
 - Platform: RC-4,1 (Code 0035) p. 25
- Implant: Ø 4,8
 - Platform: RC-4,8 (Code 0035) p. 25
- BLX
 - Implant: Ø 3,5/3,75/4/4,5
 - Platform: RB (Regular Base) (Code 0207) p. 29

TITANIUM-FIX

- b-fix
 - Implant: Ø 3,5/4
 - Platform: Regular (Code 0004) p. 22

TRE-OSS

- Simple
 - Implant: Ø 3,3/3,75/5
 - Platform: 3,75 Amarillo (Code 0040) p. 26

TRI DENTAL IMPLANTS

- TRI-Vent
 - Implant: Ø 3,75/4,1/4,7
 - Platform: 3,5 (Code 0040) p. 26

UFIT

- Gt2
 - Implant: Ø 3,5
 - Platform: Mini (Code 0004) p. 22
- Nt2
 - Implant: Ø 3,5
 - Platform: Mini (Code 0004) p. 22

VULKAN IMPLANTS

- IN-Hex
 - Implant: Ø 3,3/3,75/4,2/5
 - Platform: 3,75 (Code 0040) p. 26

YES IMPLANT

- S-SYSTEM
 - Implant: Ø 3,3/3,5
 - Platform: Narrow (Code 0030) p. 25
- Implant: Ø 4/4,5
 - Platform: Regular (Code 0030) p. 25
- Implant: Ø 5/5,5
 - Platform: Wide (Code 0030) p. 25


ZIMMER

- Screw-Vent
 - Implant: Ø 3,7/4,1
 - Platform: 3,5 (Code 0040) p. 26

COMPATIBLE with 0001

MULTI-UNIT				
	GINGIVAL HEIGHT	GINGIVAL HEIGHT	GINGIVAL HEIGHT	GINGIVAL HEIGHT
	1 mm	2 mm	3 mm	4 mm
R	42.302.001.01-2	42.302.001.02-2	42.302.001.03-2	42.302.001.04-2


DAS LIBRARIES:
MU SCANBODY DAS_MU_4.5_0001_PREMIUM
DAS_MU_8_0001_PREMIUM
DAS_MU_DS8_0001_PREMIUM
MU SCANALOG DAS_MU_SA_0001_PREMIUM
DYNAMIC μSCANBODY DAS_MU_I/IG_0001_PREMIUM

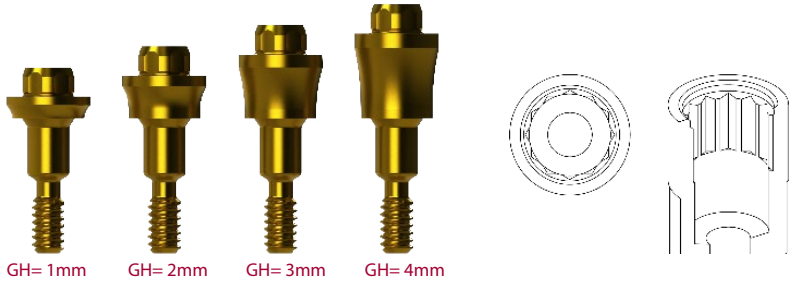

COMPATIBLE with 0002

MULTI-UNIT				
	GINGIVAL HEIGHT	GINGIVAL HEIGHT	GINGIVAL HEIGHT	GINGIVAL HEIGHT
	1 mm	2 mm	3 mm	4 mm
R	42.303.002.01-2	42.303.002.02-2	42.303.002.03-2	42.303.002.04-2


DAS LIBRARIES:
MU SCANBODY DAS_MU_4.5_0002_PREMIUM
DAS_MU_8_0002_PREMIUM
DAS_MU_DS8_0002_PREMIUM
MU SCANALOG DAS_MU_SA_0002_PREMIUM
DYNAMIC μSCANBODY DAS_MU_I/IG_0002_PREMIUM

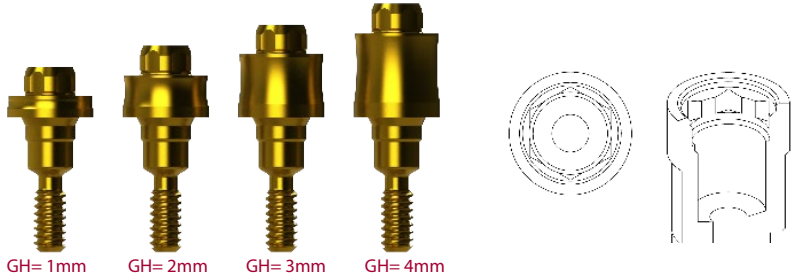

COMPATIBLE with 0004

MULTI-UNIT				
	GINGIVAL HEIGHT	GINGIVAL HEIGHT	GINGIVAL HEIGHT	GINGIVAL HEIGHT
	1 mm	2 mm	3 mm	4 mm
R	42.303.004.01-2	42.303.004.02-2	42.303.004.03-2	42.303.004.04-2


DAS LIBRARIES:
MU SCANBODY DAS_MU_4.5_0004_PREMIUM
DAS_MU_8_0004_PREMIUM
DAS_MU_DS8_0004_PREMIUM
MU SCANALOG DAS_MU_SA_0004_PREMIUM
DYNAMIC μSCANBODY DAS_MU_I/IG_0004_PREMIUM

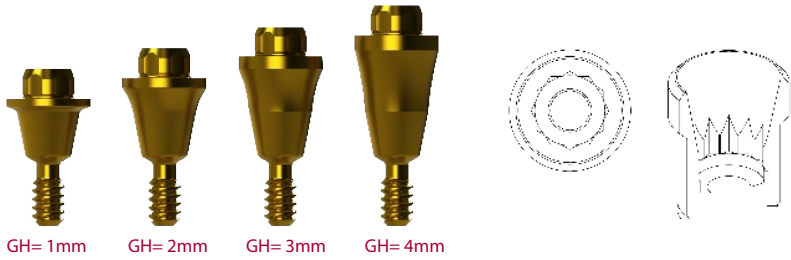

COMPATIBLE with 0006

MULTI-UNIT				
	GINGIVAL HEIGHT	GINGIVAL HEIGHT	GINGIVAL HEIGHT	GINGIVAL HEIGHT
	1 mm	2 mm	3 mm	4 mm
R	42.302.006.01-2	42.302.006.02-2	42.302.006.03-2	42.302.006.04-2


DAS LIBRARIES:
MU SCANBODY DAS_MU_4.5_0006_PREMIUM
DAS_MU_8_0006_PREMIUM
DAS_MU_DS8_0006_PREMIUM
MU SCANALOG DAS_MU_SA_0006_PREMIUM
DYNAMIC μSCANBODY DAS_MU_I/IG_0006_PREMIUM

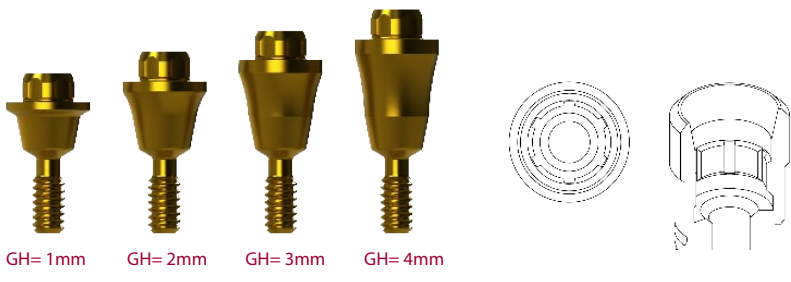

COMPATIBLE with 0007

MULTI-UNIT				
	GINGIVAL HEIGHT	GINGIVAL HEIGHT	GINGIVAL HEIGHT	GINGIVAL HEIGHT
	1 mm	2 mm	3 mm	4 mm
R	42.303.007.01-2	42.303.007.02-2	42.303.007.03-2	42.303.007.04-2


DAS LIBRARIES:
MU SCANBODY DAS_MU_4.5_0007_PREMIUM
DAS_MU_8_0007_PREMIUM
DAS_MU_DS8_0007_PREMIUM
MU SCANALOG DAS_MU_SA_0007_PREMIUM
DYNAMIC μSCANBODY DAS_MU_I/IG_0007_PREMIUM

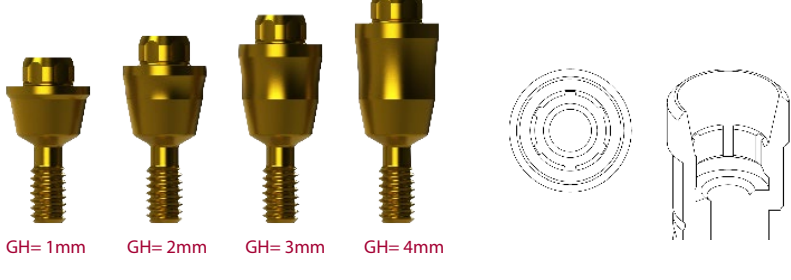

COMPATIBLE with 0014

MULTI-UNIT				
	GINGIVAL HEIGHT	GINGIVAL HEIGHT	GINGIVAL HEIGHT	GINGIVAL HEIGHT
	1 mm	2 mm	3 mm	4 mm
R	-	42.302.014.02-2	42.302.014.03-2	42.302.014.04-2


DAS LIBRARIES:
MU SCANBODY DAS_MU_4.5_0014_PREMIUM
DAS_MU_8_0014_PREMIUM
DAS_MU_DS8_0014_PREMIUM
MU SCANALOG DAS_MU_SA_0014_PREMIUM
DYNAMIC μSCANBODY DAS_MU_I/IG_0014_PREMIUM

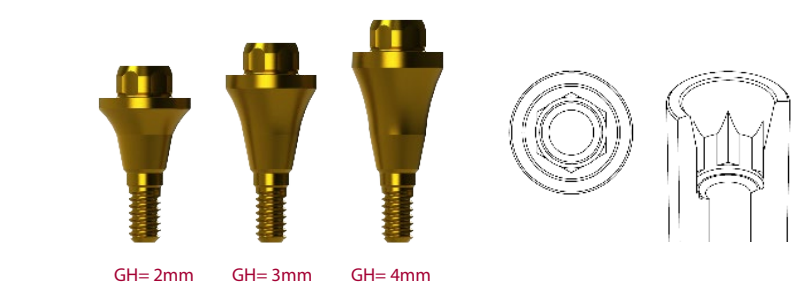

COMPATIBLE with 0015

MULTI-UNIT				
	GINGIVAL HEIGHT	GINGIVAL HEIGHT	GINGIVAL HEIGHT	GINGIVAL HEIGHT
	1 mm	2 mm	3 mm	4 mm
R	42.303.015.01-2	42.303.015.02-2	42.303.015.03-2	42.303.015.04-2


DAS LIBRARIES:
MU SCANBODY DAS_MU_4.5_0015_PREMIUM
DAS_MU_8_0015_PREMIUM
DAS_MU_DS8_0015_PREMIUM
MU SCANALOG DAS_MU_SA_0015_PREMIUM
DYNAMIC μSCANBODY DAS_MU_I/IG_0015_PREMIUM

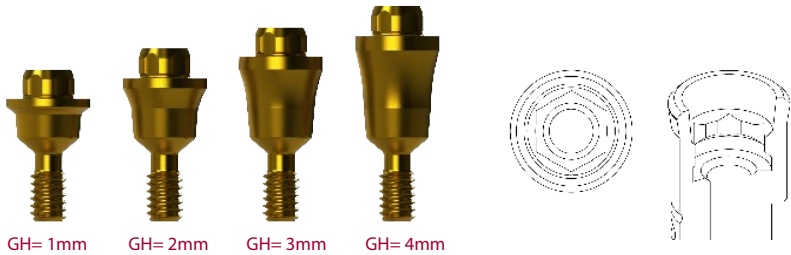

COMPATIBLE with 0021

MULTI-UNIT				
	GINGIVAL HEIGHT	GINGIVAL HEIGHT	GINGIVAL HEIGHT	GINGIVAL HEIGHT
	1 mm	2 mm	3 mm	4 mm
R	42.302.021.01-2	42.302.021.02-2	42.302.021.03-2	42.302.021.04-2


DAS LIBRARIES:
MU SCANBODY DAS_MU_4.5_0021_PREMIUM
DAS_MU_8_0021_PREMIUM
DAS_MU_DS8_0021_PREMIUM
MU SCANALOG DAS_MU_SA_0021_PREMIUM
DYNAMIC μSCANBODY DAS_MU_I/IG_0021_PREMIUM

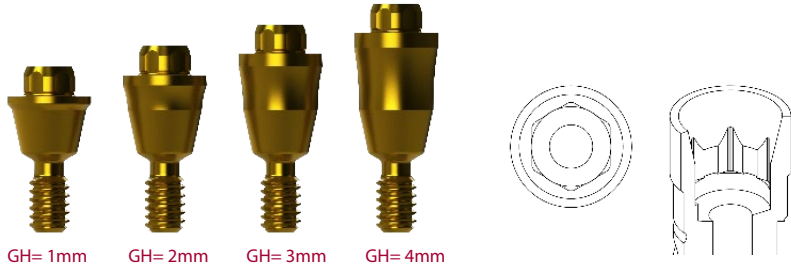

COMPATIBLE with 0022

MULTI-UNIT				
	GINGIVAL HEIGHT	GINGIVAL HEIGHT	GINGIVAL HEIGHT	GINGIVAL HEIGHT
	1 mm	2 mm	3 mm	4 mm
R	42.303.022.01-2	42.303.022.02-2	42.303.022.03-2	42.303.022.04-2


DAS LIBRARIES:
MU SCANBODY DAS_MU_4.5_0022_PREMIUM
DAS_MU_8_0022_PREMIUM
DAS_MU_DS8_0022_PREMIUM
MU SCANALOG DAS_MU_SA_0022_PREMIUM
DYNAMIC μSCANBODY DAS_MU_I/IG_0022_PREMIUM

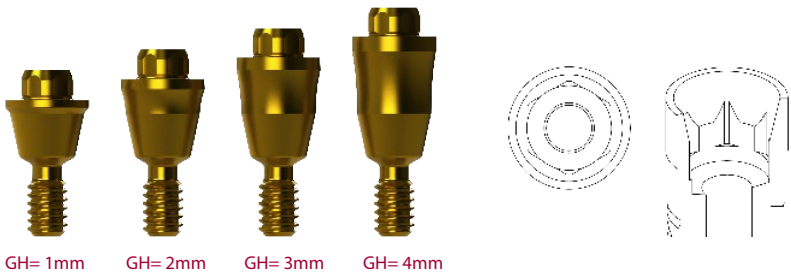

COMPATIBLE with 0029

MULTI-UNIT				
	GINGIVAL HEIGHT	GINGIVAL HEIGHT	GINGIVAL HEIGHT	GINGIVAL HEIGHT
	1 mm	2 mm	3 mm	4 mm
R	42.302.029.01-2	42.302.029.02-2	42.302.029.03-2	42.302.029.04-2


DAS LIBRARIES:
MU SCANBODY DAS_MU_4.5_0029_PREMIUM
DAS_MU_8_0029_PREMIUM
DAS_MU_DS8_0029_PREMIUM
MU SCANALOG DAS_MU_SA_0029_PREMIUM
DYNAMIC μSCANBODY DAS_MU_I/IG_0029_PREMIUM

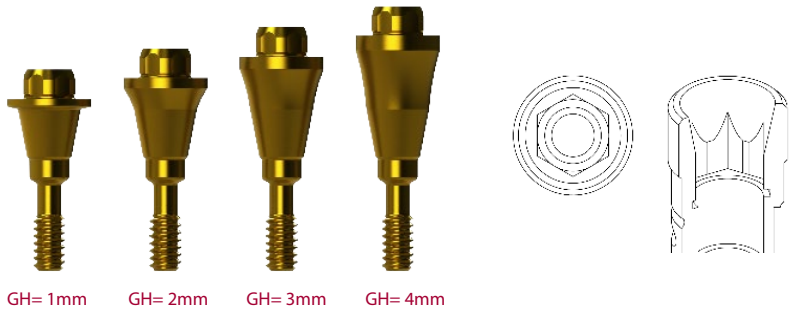

COMPATIBLE with 0030

MULTI-UNIT				
	GINGIVAL HEIGHT	GINGIVAL HEIGHT	GINGIVAL HEIGHT	GINGIVAL HEIGHT
	1 mm	2 mm	3 mm	4 mm
R	42.303.030.01-2	42.303.030.02-2	42.303.030.03-2	42.303.030.04-2


DAS LIBRARIES:
MU SCANBODY DAS_MU_4.5_0030_PREMIUM
DAS_MU_8_0030_PREMIUM
DAS_MU_DS8_0030_PREMIUM
MU SCANALOG DAS_MU_SA_0030_PREMIUM
DYNAMIC μSCANBODY DAS_MU_I/IG_0030_PREMIUM

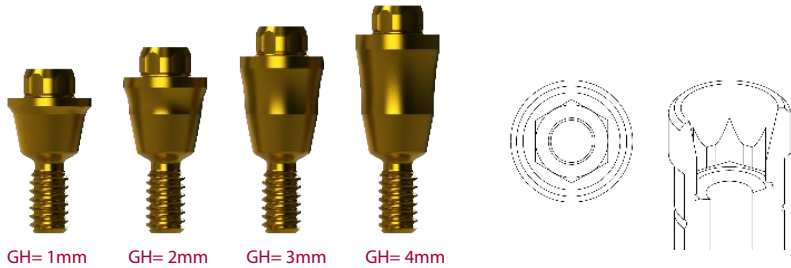

COMPATIBLE with 0033

MULTI-UNIT				
	GINGIVAL HEIGHT	GINGIVAL HEIGHT	GINGIVAL HEIGHT	GINGIVAL HEIGHT
	1 mm	2 mm	3 mm	4 mm
R	42.302.033.01-2	42.302.033.02-2	42.302.033.03-2	42.302.033.04-2


DAS LIBRARIES:
MU SCANBODY DAS_MU_4.5_0033_PREMIUM
DAS_MU_8_0033_PREMIUM
DAS_MU_DS8_0033_PREMIUM
MU SCANALOG DAS_MU_SA_0033_PREMIUM
DYNAMIC μSCANBODY DAS_MU_I/IG_0033_PREMIUM

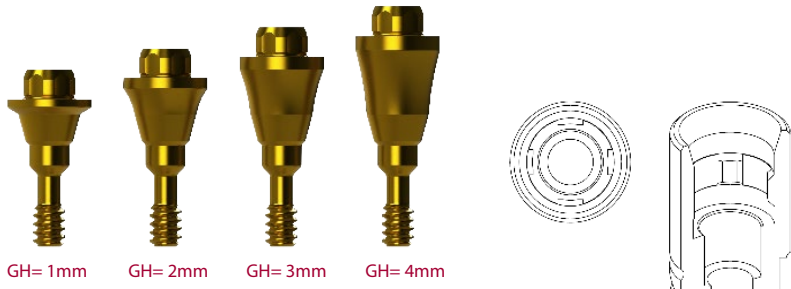

COMPATIBLE with 0035

MULTI-UNIT				
	GINGIVAL HEIGHT	GINGIVAL HEIGHT	GINGIVAL HEIGHT	GINGIVAL HEIGHT
	1 mm	2 mm	3 mm	4 mm
R	42.303.035.01-2	42.303.035.02-2	42.303.035.03-2	42.303.035.04-2


DAS LIBRARIES:
MU SCANBODY DAS_MU_4.5_0035_PREMIUM
DAS_MU_8_0035_PREMIUM
DAS_MU_DS8_0035_PREMIUM
MU SCANALOG DAS_MU_SA_0035_PREMIUM
DYNAMIC μSCANBODY DAS_MU_I/IG_0035_PREMIUM

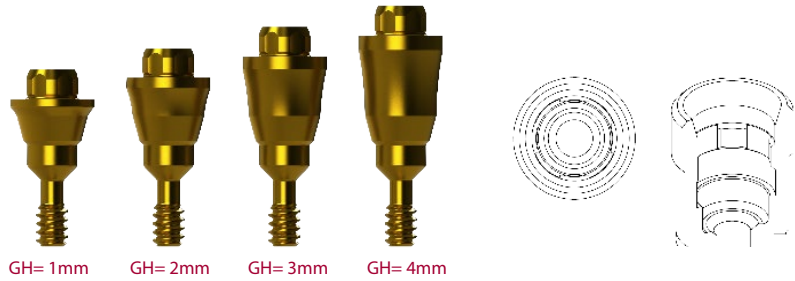

COMPATIBLE with 0040

MULTI-UNIT				
	GINGIVAL HEIGHT	GINGIVAL HEIGHT	GINGIVAL HEIGHT	GINGIVAL HEIGHT
	1 mm	2 mm	3 mm	4 mm
R	42.302.040.01-2	42.302.040.02-2	42.302.040.03-2	42.302.040.04-2


DAS LIBRARIES:
MU SCANBODY DAS_MU_4.5_0040_PREMIUM
DAS_MU_8_0040_PREMIUM
DAS_MU_DS8_0040_PREMIUM
MU SCANALOG DAS_MU_SA_0040_PREMIUM
DYNAMIC μSCANBODY DAS_MU_I/IG_0040_PREMIUM

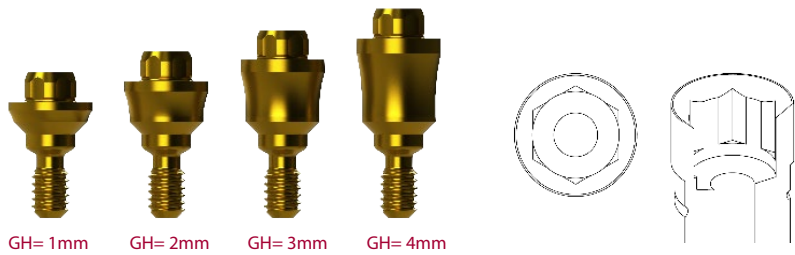

COMPATIBLE with 0075

MULTI-UNIT				
	GINGIVAL HEIGHT	GINGIVAL HEIGHT	GINGIVAL HEIGHT	GINGIVAL HEIGHT
	1 mm	2 mm	3 mm	4 mm
R	-	42.302.075.02-2	42.302.075.03-2	42.302.075.04-2


DAS LIBRARIES:
MU SCANBODY DAS_MU_4.5_0075_PREMIUM
DAS_MU_8_0075_PREMIUM
DAS_MU_DS8_0075_PREMIUM
MU SCANALOG DAS_MU_SA_0075_PREMIUM
DYNAMIC μSCANBODY DAS_MU_I/IG_0075_PREMIUM

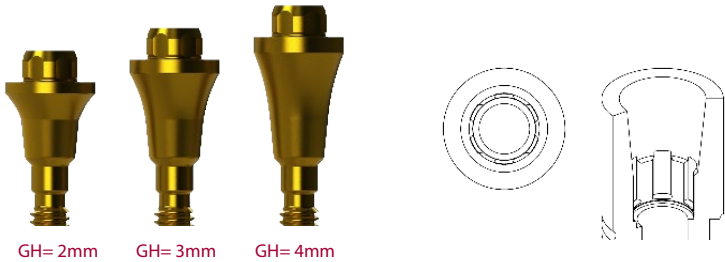

COMPATIBLE with 0091

MULTI-UNIT				
	GINGIVAL HEIGHT	GINGIVAL HEIGHT	GINGIVAL HEIGHT	GINGIVAL HEIGHT
	1 mm	2 mm	3 mm	4 mm
R	42.303.091.01-2	42.303.091.02-2	42.303.091.03-2	42.303.091.04-2

DAS LIBRARIES:
MU SCANBODY DAS_MU_4.5_0091_PREMIUM
DAS_MU_8_0091_PREMIUM
DAS_MU_DS8_0091_PREMIUM
MU SCANALOG DAS_MU_SA_0091_PREMIUM
DYNAMIC μSCANBODY DAS_MU_I/IG_0091_PREMIUM

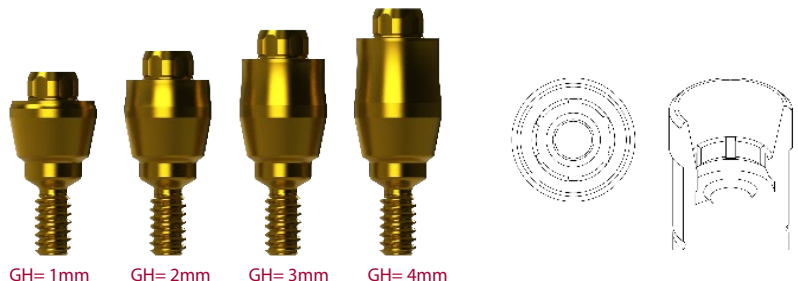

COMPATIBLE with 0102

MULTI-UNIT				
	GINGIVAL HEIGHT	GINGIVAL HEIGHT	GINGIVAL HEIGHT	GINGIVAL HEIGHT
	1 mm	2 mm	3 mm	4 mm
R	42.302.102.01-2	42.302.102.02-2	42.302.102.03-2	42.302.102.04-2


DAS LIBRARIES:
MU SCANBODY DAS_MU_4.5_0102_PREMIUM
DAS_MU_8_0102_PREMIUM
DAS_MU_DS8_0102_PREMIUM
MU SCANALOG DAS_MU_SA_0102_PREMIUM
DYNAMIC μSCANBODY DAS_MU_I/IG_0102_PREMIUM



COMPATIBLE with 0120

MULTI-UNIT				
	GINGIVAL HEIGHT	GINGIVAL HEIGHT	GINGIVAL HEIGHT	GINGIVAL HEIGHT
	1 mm	2 mm	3 mm	4 mm
R	42.303.121.01-2	42.303.121.02-2	42.303.121.03-2	42.303.121.04-2


DAS LIBRARIES:
MU SCANBODY DAS_MU_4.5_0120_PREMIUM
DAS_MU_8_0120_PREMIUM
DAS_MU_DS8_0120_PREMIUM
MU SCANALOG DAS_MU_SA_0120_PREMIUM
DYNAMIC μSCANBODY DAS_MU_I/IG_0120_PREMIUM

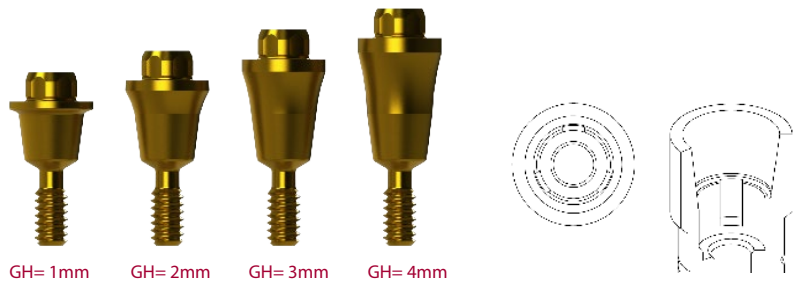

COMPATIBLE with 0121

MULTI-UNIT				
	GINGIVAL HEIGHT	GINGIVAL HEIGHT	GINGIVAL HEIGHT	GINGIVAL HEIGHT
	1 mm	2 mm	3 mm	4 mm
R	42.303.121.01-2	42.303.121.02-2	42.303.121.03-2	42.303.121.04-2


DAS LIBRARIES:
MU SCANBODY DAS_MU_4.5_0121_PREMIUM
DAS_MU_8_0121_PREMIUM
DAS_MU_DS8_0121_PREMIUM
MU SCANALOG DAS_MU_SA_0121_PREMIUM
DYNAMIC μSCANBODY DAS_MU_I/IG_0121_PREMIUM

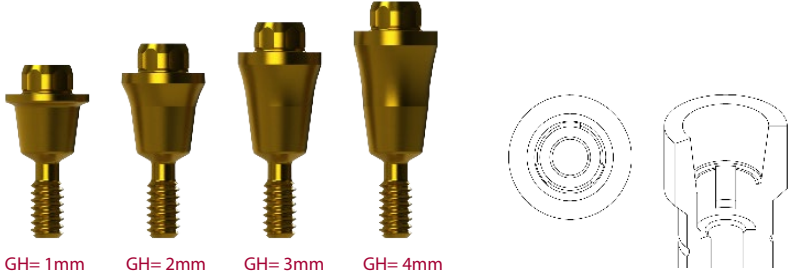

COMPATIBLE with 0125

MULTI-UNIT				
	GINGIVAL HEIGHT	GINGIVAL HEIGHT	GINGIVAL HEIGHT	GINGIVAL HEIGHT
	1 mm	2 mm	3 mm	4 mm
R	42.303.125.01-2	42.303.125.02-2	42.303.125.03-2	42.303.125.04-2

DAS LIBRARIES:
MU SCANBODY DAS_MU_4.5_0125_PREMIUM
DAS_MU_8_0125_PREMIUM
DAS_MU_DS8_0125_PREMIUM
MU SCANALOG DAS_MU_SA_0125_PREMIUM
DYNAMIC μSCANBODY DAS_MU_I/IG_0125_PREMIUM

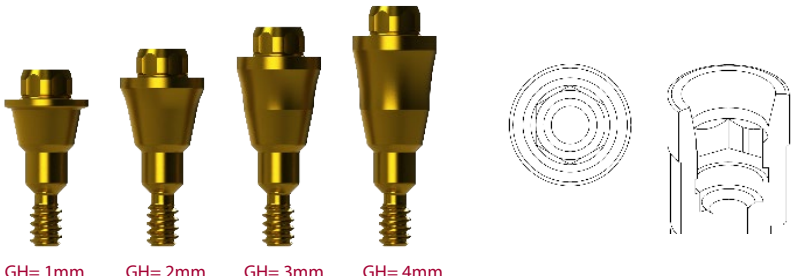

COMPATIBLE with 0186

MULTI-UNIT				
	GINGIVAL HEIGHT	GINGIVAL HEIGHT	GINGIVAL HEIGHT	GINGIVAL HEIGHT
	1 mm	2 mm	3 mm	4 mm
R	42.303.186.01-2	42.303.186.02-2	42.303.186.03-2	42.303.186.04-2

DAS LIBRARIES:
MU SCANBODY DAS_MU_4.5_0186_PREMIUM
DAS_MU_8_0186_PREMIUM
DAS_MU_DS8_0186_PREMIUM
MU SCANALOG DAS_MU_SA_0186_PREMIUM
DYNAMIC μSCANBODY DAS_MU_I/IG_0186_PREMIUM



COMPATIBLE with 0169

MULTI-UNIT				
	GINGIVAL HEIGHT	GINGIVAL HEIGHT	GINGIVAL HEIGHT	GINGIVAL HEIGHT
	1 mm	2 mm	3 mm	4 mm
R	42.302.169.01-2	42.302.169.02-2	42.302.169.03-2	42.302.169.04-2

DAS LIBRARIES:
MU SCANBODY DAS_MU_0169_PREMIUM
DAS_MU_8_0169_PREMIUM
DAS_MU_DS8_0169_PREMIUM
MU SCANALOG DAS_MU_SA_0169_PREMIUM
DYNAMIC μSCANBODY DAS_MU_I/IG_0169_PREMIUM



COMPATIBLE with 0207

MULTI-UNIT				
	GINGIVAL HEIGHT	GINGIVAL HEIGHT	GINGIVAL HEIGHT	GINGIVAL HEIGHT
	1 mm	2 mm	3 mm	4 mm
R	42.302.207.01-2	42.302.207.02-2	42.302.207.03-2	42.302.207.04-2

DAS LIBRARIES:
MU SCANBODY DAS_MU_0207_PREMIUM
DAS_MU_8_0207_PREMIUM
DAS_MU_DS8_0207_PREMIUM
MU SCANALOG DAS_MU_SA_0207_PREMIUM
DYNAMIC μSCANBODY DAS_MU_I/IG_0207_PREMIUM



TALLADIUM GUARANTEE

TERMS AND CONDITIONS

These guarantee terms and conditions ("T&C") cover the entire range of Talladium products ("Products"), manufactured by TALLADIUM ESPAÑA S.L. and distributed by Geoda Medical S.L. or official dealers. The guarantee described in these T&C is exclusively in benefit of the clinician ("Clinician") and of the dental technician ("Technician") and not for the benefit of third parties or institutions, including patients.

GUARANTEE PERIOD

TALLADIUM ESPAÑA S.L. offers a lifelong guarantee for its entire range of products starting from the date of issue of the invoice.

GUARANTEE SCOPE

Subject to the limitations and exceptions described in these T&C, TALLADIUM ESPAÑA S.L. will offer the following benefits:

QUALITY: If there are defects in the materials or in the manufacturing of the Product, TALLADIUM ESPAÑA S.L. will replace the Product with no additional cost.

SAFETY: If, having complied with all the product indications, the prosthesis should have to be made again, due to a fault in the Dynamic Abutment® or Dynamic Titanium Base® system, TALLADIUM ESPAÑA S.L. will replace the abutments and screws necessary to remake the prosthesis, as well as the costs derived from its manufacturing.

In case of having used our products and having complied with all the product indications, the implants suffer any damage, TALLADIUM ESPAÑA S.L. will pay the cost of the implants. This coverage will only be valid during the first 6 months after the collocation of the prosthesis which includes our products.

CLAIM REQUIREMENTS AND PROCEDURE

To receive the benefits indicated in these T&C, the treating Clinician must satisfy the following requirements:

- a) The claim must be notified to TALLADIUM ESPAÑA S.L. within (30) days since the date the claimed defect was detected.
- b) This requires that the Clinician or Technician must contact the customer service department by telephone or by e-mail to make the claim.
- c) A claim form will be completed, which, together with a document or report which justifies the faulty Product and the faulty Product itself, will be sent by the customer to TALLADIUM ESPAÑA S.L. offices, within the previously indicated period.
- d) Clinicians or Technicians presenting a claim in agreement with these T&C must be up to date in any payments owing to TALLADIUM ESPAÑA S.L. or to any of its subsidiaries, at the time when the claim form is presented.
- e) All the use procedures of our Products must be carried out in agreement with the instructions of TALLADIUM ESPAÑA S.L. as well as in accordance with commonly accepted dentistry practices.
- f) The expenses derived from this procedure will be assumed by the customer. The return shipping costs will be assumed by TALLADIUM ESPAÑA S.L. in all those cases covered by these T&C.

Regardless of the guarantee rights, claims should be notified as soon as possible in order to comply with regulatory requirements.

GENERAL LIMITATIONS OF THIS GUARANTEE

With the exception of the guarantee described in these T&C, neither TALLADIUM ESPAÑA S.L. nor its representatives, nor third parties manufacturing or distributing the Products, represent or offer a guarantee, agreement or any other express or implicit, oral or written, commitment, with respect to the Products (without limitation), including guarantees involved in the marketing, durability or suitability for individual uses or purposes.

In addition and within the maximum extent permitted by the relative law, TALLADIUM ESPAÑA S.L. rejects (on its own behalf, and on behalf of its representatives and third parties that manufacture or distribute Products) any responsibility with respect to any direct or indirect damage caused, which may result from or be a consequence of the design, composition of the dental prosthesis into which the Products are integrated.

GUARANTEE EXCLUSIONS

TALLADIUM ESPAÑA S.L. limits this guarantee to:

- Transformed abutments that form part of the dental prosthesis. But not the screws used to anchor them.
- Clinical screws that have been in the mouth for more than 2 years.

AMENDMENT OR SUSPENSION OF THE GUARANTEE

TALLADIUM ESPAÑA S.L. reserves the right to amend or withdraw these T&C at any time and without prior notification. Any modification or suspension shall not affect products already placed in patients.

Ed.2019-01



DYNAMIC ABUTMENT® SOLUTIONS

VIRGINIA WOOLF, 17 · 25005 · LLEIDA (SPAIN)

INTERNATIONAL +34 873 450 709

das@dynamicabutment.com

SPAIN +34 973 289 580

spain@dynamicabutment.com

Dynamic Abutment® Solutions is a trademark of the manufacturer Talladium España S.L

The information included in this catalogue is exclusively addressed to professionals in dental sector.

All commercial trademarks mentioned herein are fully registered by their respective companies, and the images that appear are just to provide an orientation.

Talladium reserves the right to modify this document without previous notice.

Talladium is not responsible for the inadequate execution of these products if the warning indications corresponding to every reference are not contemplated.

All of the products listed in this brochure are marking in accordance with CE legislation. Some of the products are not authorized for sale and distribution or do not have a sales license in some countries according to other legislations (FDA, MDSAP, etc.) Please ask for information: das@dynamicabutment.com

 *Marking in accordance with CE legislation and applicable sanitary regulations*



Visit our Online Store to find all our products and compatibilities :

www.dynamicabutmentstore.com

Ed.2021-04



DYNAMIC ABUTMENT® SOLUTIONS

Virginia Woolf, 17

25005 · Lleida · Spain

www.dynamicabutment.com

INTERNATIONAL +34 873 450 709

das@dynamicabutment.com

SPAIN +34 973 289 580

spain@dynamicabutment.com

DYNAMIC ABUTMENT® SOLUTIONS